Conndevs

Awesome connected devices for everyday use

 Smart Body Composition Scale

User Manual

ENGLISH	p.1
DEUTSCH	p.2
FRANÇAIS	p.4
ITALIANO	p.5
ESPAÑOL	p.6

ENGLISH

Thank you for choosing our Smart Body Composition Scale!

BIA TECHNOLOGY

This model utilizes the BIA (Bio-electrical Impedance Analysis) technology. This cutting-edge technology makes the measurements more accurate by taking advantage of the conductive composition of muscle and water contents within the body. The scale sends harmless signals through the body to measure the overall body impedance, which is increased by fat composition (not conductor) and lowered by the muscles and water in the body (conductors)

PRODUCT PARAMETERS


- 1.Product Size: 280*280*23mm
- 2.Product Weight: 1.2kg/ 2.64lb
- 3.Unit conversion: kg/ lb
- 4.Range: 5kg/ 11lb – 180kg/ 397lb
- 5.Division value: 0.1kg/ 0.2lb
- 6.Age: 10 years old or older
- 7.Height range: 100 ~ 220 cm
- 8.Working temperature/ humidity: -10°C ~ 50°C/ 20% ~ 90%
- 9.Battery installation: 3X1.5V AAA
- 10.WiFi Band: 2.4 GHz only

PRECAUTIONS & SECURITY

- 1.Pregnant women and the users wearing pacemakers or other dangerous devices are not recommended to use this product.
- 2.This scale should not be used to diagnose or treat any medical condition. You should consult a primary care doctor or physician before undertaking your diet or exercise program.
- 3.Do not overload the scale (maximum 180kg/ 397lb)
- 4.Do not immerse the scale on water
- 5.Please use on a flat hard floor, do not use on slippery floors.
- 6.Do not stand on the edge of the scale or jump on it, in order to avoid falling.
- 7.Children under the age of 10 should not touch the product to avoid injury.
- 8.Keep your feet dry and clean when stepping on the scale.
- 9.In order to obtain accurate results, please measure at the same time every day.
- 10.Do not test after strenuous exercise.
- 11.Do not disassemble this product yourself.
- 12.Do not scrub this product with water or volatile oil.
- 13.Please put the batteries in the correct position. If not used for a long time, please take out the battery.
- 14.Test results are for reference only, and please do not use as a basis for your diet or fitness program.

WIFI CONNECTION & HEALTHU + APP

Please download the HealthU + app from the Apple App Store or Google Play Store. Then follow the connection instructions. At a certain point, you'll have to select the WiFi Scale (instead of Bluetooth). Please follow app instruction for connecting the scale to WiFi and the app. HealthU + app was designed to be very intuitive and easy to use. You can set the language from the Settings Menu. If you need assistance, you can also check the FAQs on the Settings Menu

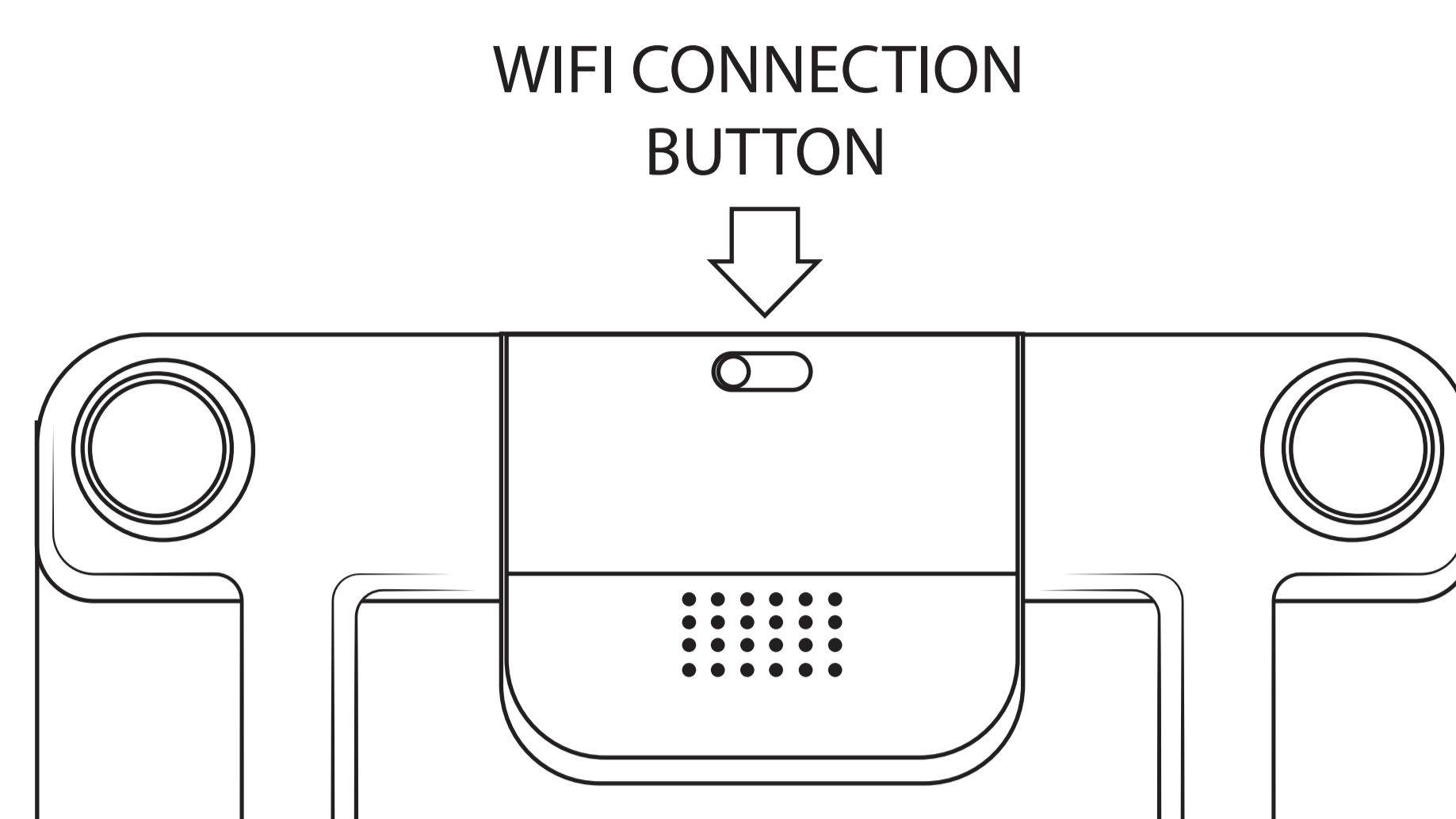FAQ & CUSTOMER SUPPORT

Please refer to the FAQs on the Settings Menu of the app. In case you need further assistance, please contact us at support@conndevs.com. We will be happy to answer all your questions and comments!

WHAT'S INCLUDED

Smart WiFi Scale / Batteries: 3X1.5V AAA / User Manual

GUARANTEE

The scale has passed quality controls and has been tested for accuracy prior to leaving the factory. It has a guarantee against any defect in either workmanship or materials for a period of one year.

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